


Cuspidor Valve for Installation in the Suction Pipe (below the Cuspidor)

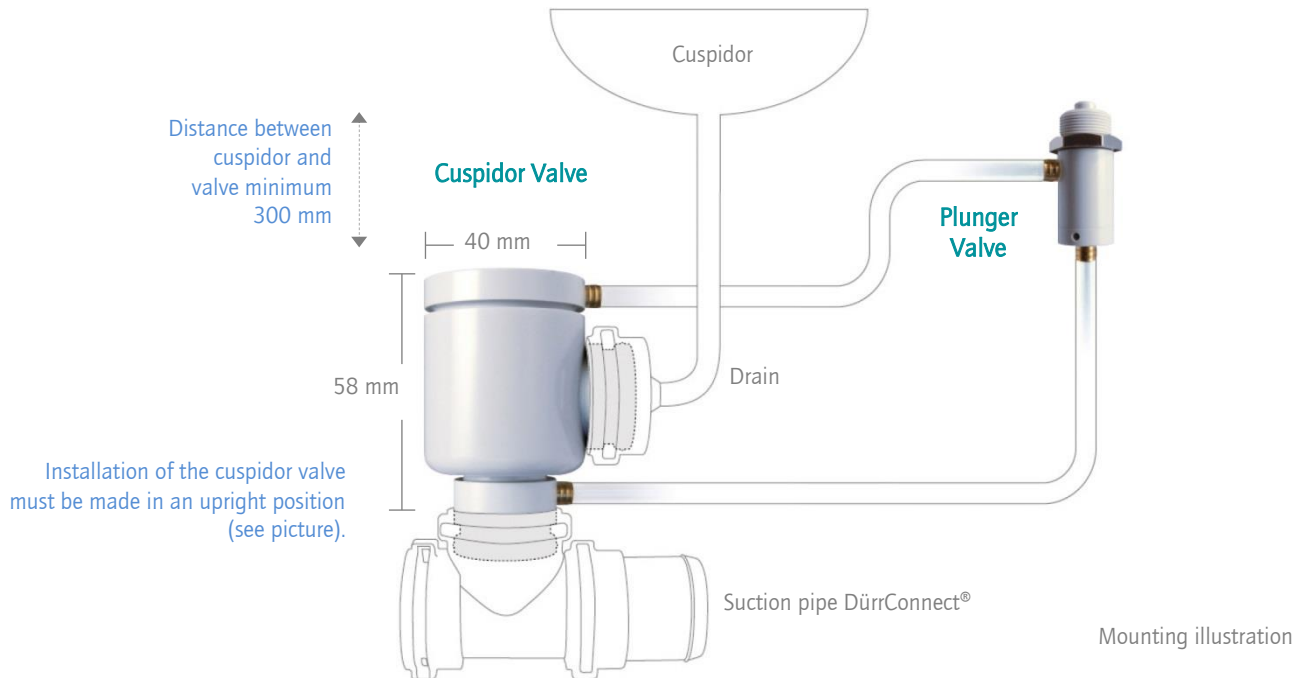

! Please note that

- installation of the cuspidor valve must be made in an upright position (see picture).
- distance between cuspidor and valve must be at least 300 mm.
- suction performance has to be between 300 and 350 l/min.
- the suction machine must be switched on by removing the suction hose from the hose manifolds.

One contact only!

We have our own development department, manufacture our own tools and moulds and produce over 90 % of our products ourselves. Therefore, we are in a position to respond to any new product requirements quickly and can ensure highest quality standards. Medtronic is certified according to DIN EN ISO 13485.