


Water Bottle System Clean Water

for mounting ON unit (via angle) or IN unit (mounting in unit from below)

for mounting ON type with mounting bracket stainless steel or white
 with thread for S-Stream bottles



with standard thread for e. g. DCI® bottles

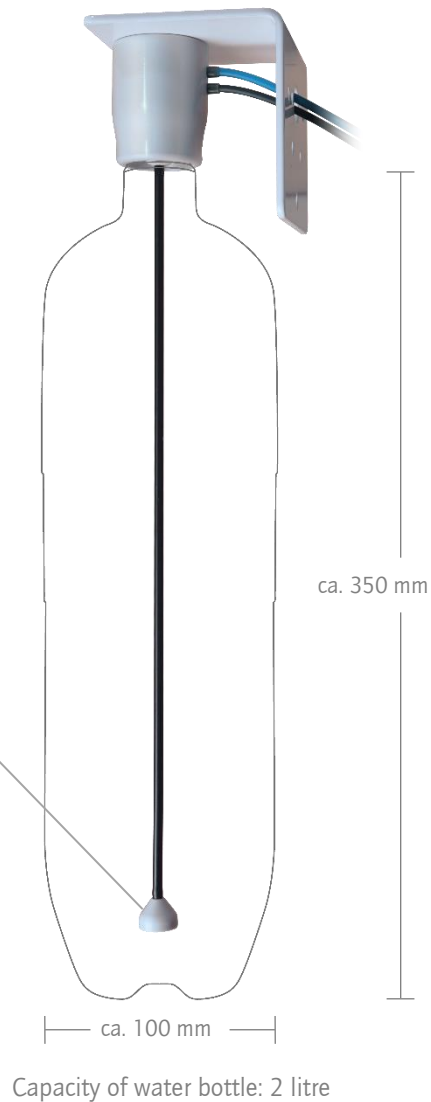

! For hygienic safety you can rely on Medtronic's Clean Water System

- Can be adapted to any treatment unit
- Use of clean water ensures technical safety

For further details please see reverse!

Water Bottle System

For Hygienic Safety You Can Rely on Medtronic's Clean Water System

- The system can be adapted to any treatment unit
 - Mounting on treatment unit by using a mounting bracket
 - Water bottle connected to compressed air (regulated by an automatic pressure regulator)
 - Automatic pressure release when exchanging the water bottle
- The use of clean water ensures technical safety
 - Medtronic's water bearing hoses are approved for drinking water use and are biocompatible.
 - Our hoses are tested in compliance with DVGW worksheet no. W 270*. (Reduction of micro-biological growth such as bio-slime, fungi and bacteria in drinking water applications).


* DVGW = German Association of Gas and Water

- Stainless steel version for mounting ON type with or without mounting bracket




Mounting IN type

with thread for S-Stream bottles

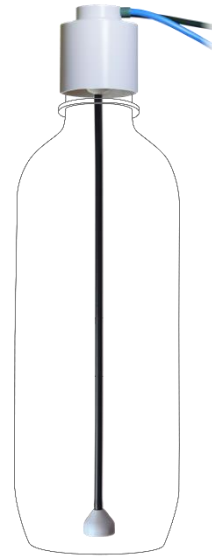

with standard thread for e.g. DCI® bottles

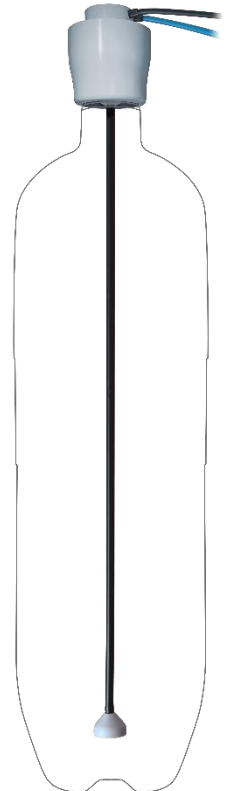

! The Medtronic Water Bottle Systems are also available with a higher water flow rate for cup filling! Please contact us for further information!



medizinisch-elektronische Geräte-Gesellschaft m.b.H.

Riedelbach | Technologiepark
In den Ensterwiesen 7
61276 Weilrod | Germany

Tel +49 (0) 6083 95989-0
Fax +49 (0) 6083 95989-11
info@medtronic-dental.de
www.medtronic-dental.de

One contact only!

We have our own development department, manufacture our own tools and moulds and produce over 90 % of our products ourselves. Therefore, we are in a position to respond to any new product requirements quickly and can ensure highest quality standards. Medtronic is certified according to DIN EN ISO 13485.